



AURORA[®]
special stroller



Aurora® is an advanced special needs stroller that combines a new generation design with a modern look. It is designed for the youngest patients with limited locomotor abilities. It holds the most essential features of the special needs stroller, not resembling one at the same time. The proven technical solutions used in the stroller make it safe, stable and very comfortable. The biggest advantage of the **Aurora®** stroller is its convenient multifunctional reversible seat. It allows the child to interact with the caregiver or the environment, depending on the seat position. A tilt-in-space function allows for smooth adjustment of the inclination angle of the entire seat. The adjustable cushions and side supports provide excellent lateral stabilization. In addition, there is the possibility of changing the height of the headrest, the angle of the backrest and the footrest. The **Aurora®** stroller has been equipped with adjustable shock absorbers. There are two modes to choose from: Soft and Hard. The front swivel wheels with blockade make it easier to maneuver the stroller in tight and narrow spaces. As standard, the stroller is equipped with practical accessories: an adjustable folding canopy, leg cover, capacious under seat storage basket and bag. Numerous reflective elements placed on the stroller provide additional safety. The stroller is available in one universal size. A wide range of adjustments make it possible to adapt the stroller to the individual needs of smaller or larger children.



IT IS RECOMMENDED FOR CHILDREN IN THE FOLLOWING CASES:

- Infantile cerebral palsy – different forms.
- Meningomyelocele.
- Brain and spinal cord injuries.
- Muscle diseases connected with paresis and paralysis.
- Genetic disorders connected with paresis.
- Other diseases connected with paresis, paralysis and injuries regarding locomotor system.
- After traumatic brain injury, spine with spinal cord and lower limbs injury.



STANDARD EQUIPMENT:

Foldable frame | Hard/soft suspension system | Detachable and reversible seat unit | Fully reclining backrest converting the stroller into a lounger (couch function) | Tilt in-space adjustment | Front PU swivel wheels with direction blockade - quick release | Rear PU wheels with central brake - quick release | Adjustable footrest with supporting belt | Adjustable push handle | Under seat storage basket | Removable abduction block | Adjustable lateral supports | Seat narrowing pads | 5 points safety belts with a Velcro fastening | Adjustable headrest | Folding canopy | Leg cover | Front rail | Rain cover | Bag



ADVANTAGES:



Stable structure

The Aurora® stroller provides proper stabilization and the user safety.



Reversible seat

There is the possibility of putting the seat in the forward or rearward facing position, depending on the needs of the child.



Light weight

The Aurora® stroller is made of light materials, which is an additional benefit.



Easy-to-use

The Aurora® stroller has intuitive adjustments and its assembly is simple.



24 months warranty

We offer 24 months warranty on the product. We also provide after sales service.



Tilt-In-Space

It is possible to change the tilt angle of the seat.



The tilt-in-space function can be easily achieved by the adjustment of the frame



2 available settings of the shock absorbers: Hard and Soft



The possibility of putting the seat in the forward or rearward facing position



Adjustable backrest

COLORS

WATERPROOF UPHOLSTERY HYDRO™



WATERPROOF



EASY TO CLEAN



DURABLE



ECOLOGICAL



FIRE RETARDANT



Melange




Vulcano



Denim

* Optional equipment are custom-made. lead time ca. 30 days

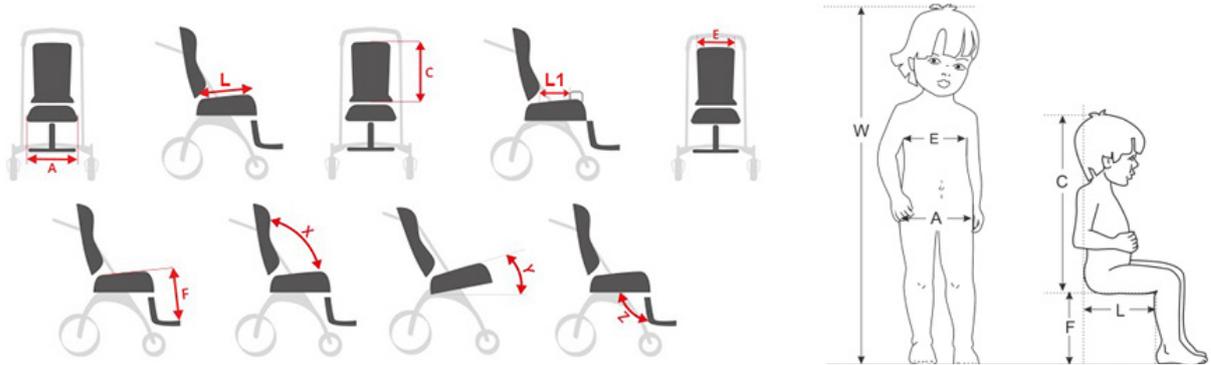

Symbol	Measurement	Unit	Size 1
W	Indicative user height***	[cm]	110
A	Seat width	[cm]	18*/27
L	Seat depth	[cm]	30
C	Headrest height	[cm]	43 - 51
E	Backrest width	[cm]	20**/28
F	Footrests length	[cm]	23, 26, 28
	Backrest height	[cm]	55
L1	Seat depth to the abduction block	[cm]	11 - 22
	Footrests length	[cm]	18,5
X	Backrest tilt angle	[°]	110 - 180
Y	Seat tilt angle	[°]	0 - 30
Z	Footrests tilt angle	[°]	90 - 180
	Maximum user weight	[kg]	40

Note: Measurement tolerance range is +/- 1 cm / ½ in.

Note: Approximate values, depending on the product configuration.

* Measurement starts from narrowing pad.

** Measurement starts from lateral support.

*** Height cannot be main factor for selection of product size.

Product dimensions & weight		Size 1
Width	[cm]	67
Length	[cm]	118
Height	[cm]	96
Width after folding	[cm]	62
Length after folding	[cm]	95
Height after folding	[cm]	47
Weight	[kg]	19

Note: Measurement tolerance range is +/- 1 cm / ½ in.

Note: Approximate values, depending on the product configuration.



Safety



ARO_001 LED lights
Front and rear LED lights for increased safety.



Positioning



ARO_160 Octagon™ Headrest
Equipped with an internal aluminum 'fin' covered with soft foam. Effectively stabilizes the patient's head through greater precision of the shape adjustment. The headrest cover is made of elastic, durable and dirt-repellent fabric and is certified Oeko-Tex 100.



ARO_161 Sandals stabilizing feet and ankles
They protect and stabilize the ankle joint. They are made of elastic material, recommended for users with increased muscle tone. They can be used with or without shoes.



ARO_162 Dynamic ankle stabilizers (FP-22)
Dynamic FlexPoint® ankle stabilizers provide great support when stabilizing the user's feet. Incorrect positioning of the feet, despite appropriate stabilization of the pelvis and chest, makes it impossible to keep the correct sitting position.



Comfort



ARO_415 Rear wheels cover
Use it on the dirt tires when you want to take the products indoor or transport in the trunk.



ARO_429 Front wheels cover
Comfortable pad, protecting the feet from direct contact with the footplate



ARO_442 Front wheels cover
Use it on the dirt tires when you want to take the products indoor or transport in the trunk.



ARO_447 Winter muffs
Comfortable and warm muffs protect your hands against a chill during walks.



ARO_448 Lap tray



Storage



ARO_509 Organizer
Handy carrying bag on the handle of the trolley. Specially designed compartments make it easy to keep order by storing the most needed items at hand. The Organizer has two useful holes for a cup of coffee or a bottle of water.