

KUKINI™

Special needs stroller



Class I Medical Device

**This is a medical device.
Use it in accordance with the user manual or the label.**

Description

To meet the expectations of our customers, the Akces-MED Ltd. team created a new special needs stroller **KUKINI™**. It's really light so it is perfect for both long walks and running trainings. Aluminium frame of the stroller with reflective elements makes the stroller safe and comfortable for your child while exploring the city or forest areas. **KUKINI™** is equipped with a rear shock absorbing system and 16-inch inflatable spoked wheels, which make riding on uneven terrain a pleasure for both a child and caregiver.

The child's safety and correct sitting position are guaranteed by a 5-point safety belts with a magnetic fastener. The comfortable head support is provided by the anatomical **Octagon™** headrest, which can be mounted at different heights. An adjustable backrest allows you to set the optimal position for your child while riding. The wheel quick release system together with the compact size of the folded stroller make it easy to store or transport. Durable and comfortable upholstery has a practical pocket to hold your phone, wallet or keys.

In order to protect the child from various weather conditions, the special needs stroller **KUKINI™** can be equipped with an adjustable folding canopy to provide additional protection from wind, rain or sun.

Indications

- Cerebral palsy – various types and forms.
- Brain and spinal cord injuries with lower limbs paralysis and weakness, balance and motor coordination disorders with different etiology.
- Demyelinating disease (e.g. multiple sclerosis – SM).
- Meningocele.
- Muscles diseases with weakness and paralysis (dystrophy, myasthenia gravis, myositis).
- Genetic disorders.
- Degeneration and neurological diseases.
- Other condition connected with muscles weakness, paralysis and musculoskeletal disorders.
- After head and spinal cord injuries and spine trauma with paralysis and paresis.



Special features



Stable structure

The device ensures adequate stability and safety of user.



Easy to use

Intuitive way of adjustment and easy assembly.



Lightweight construction

The product is made of light metal alloys which positively affects the low weight of the device.



Reflective elements

The product has reflective elements to ensure your safety.

24⁺

24-month warranty

We offer 24-month warranty on the product. We also provide after sales service.

Colors

Waterproof upholstery **HYDRO™**



Gray

You have full control.

The disc brake, operated by a central handle located on the comfortable, ergonomic push handle, allows you to control your speed. Invaluable when climbing hills.

A wrist safety strap prevents accidental rolling away of the stroller with a child. Rear wheel central brake ensures safety when the stroller is not in motion, for example when putting the child in or taking them out of the stroller.



We focused on every detail.

Cleverly designed, lightweight aluminium frame guarantees safety and comfort. Rear shock absorbing system and 16-inch inflatable spoked wheels make riding a pleasure, even on rough terrain. Easy folding system and compact size when folded.

One-of-a-kind upholstery

Durable waterproof upholstery and a comfortable seat cushion with VISmemo™ memory foam provide your child with great comfort while riding. The upholstery is equipped with reflective elements to improve visibility on the road.



A born champion.

A special needs stroller KUKINI™ is a perfect solution for people who look for comfort, safety and stability for their children, but don't want to give up an active way of spending time walking or running.



Special needs strollers



Standard equipment

Ergonomic push handle

Backrest with adjustable tilt angle and a storage pocket

5-point safety belts with a magnetic fastener

Seat cushion with VISmemo™ foam

Foldable aluminum frame equipped with a front wheel disc brake and rear wheel central brake

Removable inflatable rear wheels

Adjustable footrest

Removable inflatable front wheel

+ Voyager backpack



Important measurements

Sizes	Seat width	Seat depth	Backrest width	Backrest height	Backrest tilt angle	Footrest length
Size 1	37 cm	37 cm	36 cm	57 cm	100 - 120°	27, 40 cm

Size	Indicative user height*	Max user weight
Size 1	100 - 135 cm	40 kg

Weight and dimensions

Size	Width	Width after folding**	Length	Length after folding**	Height	Height after folding**	Weight
Size 1	67 cm	54 cm	124 - 143 cm	97 cm	105 cm	38 cm	14 kg

Scan the QR code to see the full table.



*Height cannot be the main parameter while choosing the product size.
 **Folded dimensions are for a stroller with removed wheels and an unfastened footrest


Note: Measurement tolerance range is +/- 1 cm / ½ in.
 Note: Approximate values, depending on the product configuration.

Accessories



Folding canopy

Protects against both sun and rain.

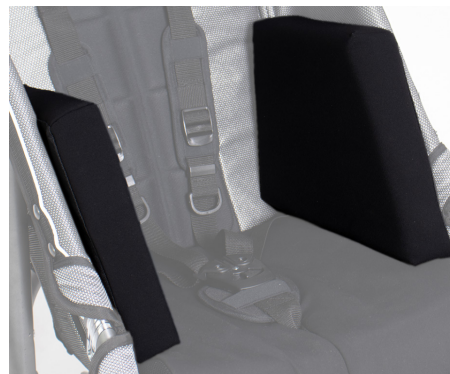

Narrowing pads 6/10cm

Made of elastic foam trimmed with soft material. They are attached to the seat upholstery from both sides. Thanks to them you can narrow the seat by 6/10 cm.



Mosquito net

Excellent protection against insects.



Bicycle tow bar with a flag

Specjalny adapter zmienia wózek KUNINI™ w przyczepkę rowerową. Wystarczy zamocować go na tylnym kole roweru i delektować się podróżami wraz z dzieckiem siedzącym bezpiecznie i wygodnie w KUKINI™.



Winter footmuff with muffs

Warmly lined footmuff protects the child from cold during winter walks and gives the user the freedom of movement. The zipper allows you to quickly unfasten the footmuff, making it easier to place the user in it.



Scan the QR code, to see all the accessories.