

# HARMONY™

Special needs stroller



**MD** Class I Medical Device

**This is a medical device.  
Use it in accordance with the user manual or the label.**

## Description

---

**HARMONY™** is a special needs stroller designed for children with musculoskeletal dysfunctions. It has the most important features of a therapeutic wheelchair, while not resembling it at all. Thanks to its proven technical solutions, it is safe, stable and at the same time very comfortable.

**HARMONY™** provides a high degree of stability and fit. The length-adjustable backrest, seat and single-column footrest allow the stroller to be properly adjusted to the child's height. The head stabilizer, along with torso and hip pelvises, support the little user's body by providing lateral stabilization. The comfortable multifunctional seat has a front- to-back option. This gives the child the option of sitting facing forward while observing the surroundings or backward while maintaining eye contact with the caregiver. The bucket function allows you to smoothly adjust the tilt angle of the entire seat.

**HARMONY™** stroller is equipped with a cushioned suspension system that is adjustable. There are two modes to choose from: City for riding on a flat, level road and Comfort for riding on difficult, uneven terrain. The front swivel wheels with lockable steering further help maneuver the stroller in tight spaces.

Available in 1 universal size. It is designed for the youngest users up to about 120 cm in height.

## Indications

---

- Infantile cerebral palsy – different forms.
- Meningomyelocele.
- Brain and spinal cord injuries.
- Muscle diseases connected with paresis and paralysis.
- Genetic disorders connected with paresis.
- Other diseases connected with paresis, paralysis and injuries regarding locomotor system.
- After traumatic brain injury, spine with spinal cord and lower limbs injury.



## Special features



### **Stable structure**

The device ensures adequate stability and safety of user.



### **Easy to use**

Intuitive way of adjustment and easy assembly.



### **Lightweight construction**

The product is made of light metal alloys which positively affects the low weight of the device.



### **Reflective elements**

The product has reflective elements to ensure your safety.



### **Tilt-In-Space**

The device is able to change hip angle.



### **Movable seat**

You can turn the seat easily so that the child may be facing towards caregiver or movement direction.

24<sup>+</sup>

### **24-month warranty**

We offer 24-month warranty on the product. We also provide after sales service.

## Colors

---

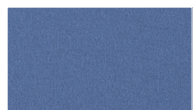
Waterproof upholstery **HYDRO™**



**MELANGE**



**VULCANO**



**COBALT**



**FLAMING**



## Movable seat

---

The seat can be put in the forward or rearward position.





## Stabilization

---

A proper stabilization of the user is provided by a wide range of standard equipment.



## Tilt-In-Space

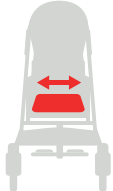







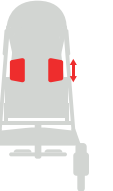





---

The tilt-in-space function can be easily achieved by the adjustment of the frame.

# Standard equipment



# Important measurements

Sizes	 Seat width (without pelvic support)	 Seat depth	 Seat height	 Backrest height	 Armrests height	 Spacing between pelvic support	 Spacing between chest support	 Chest support height	 Spacing between head support	 Footrest height	 Footplate length	 Footplate tilt angle
Size 1	34 cm	30 - 35 cm	57 cm	54 - 60 cm	12 - 17 cm	17 - 30 cm	18 - 24 cm	26 ÷ 29 MIN backrest height; 32 ÷ 35 MAX	12 - 20 cm	13 - 29 cm	18 cm	-5 - 10°
Size	 Indicative user height*		 Max user weight									
Size 1	100 - 120 cm		40 kg									

# Weight and dimensions

Size	Width	Width after folding	Length	Length after folding	Height	Height after folding	Weight
Size 1	61 cm	60 cm	98 cm	95 cm	73 - 118 cm	47 cm	21 kg

Scan the QR code to see the full table.



\*Height cannot be the main parameter while choosing the product size.

Note: Measurement tolerance range is +/- 1 cm / ½ in.  
Note: Approximate values, depending on the product configuration.

# Accessories



## 4 points safety vest

Acts as a trunk stabilizer. Wide range of adjustment helps to fit it to the user's needs.



## Dynamic ankle stabilizers (FP-22)

ABS thermoformable tray, perfect for recreation and during rehabilitation. It is also ideal for work and having meals.



## PA-VM-15 cushion

VISmemo™ memory foam used in the cushion increases the seating comfort. A bolster at the front of the cushion helps maintain the position of the pelvis, reducing the possibility of sliding off the seat.



## Sandals stabilizing feet and ankles

They protect and stabilize the ankle joint. They are made of elastic material, recommended for users with increased muscle tone.



## Tray

ABS thermoformable tray, perfect for recreation and during rehabilitation. It is also ideal for work and having meals.



Scan the QR code, to see all the accessories.