

ZEBRA INVENTO™ NEW VERSION

Positioning chair



**This is a medical device.
Use it in accordance with the user manual or the label.**

Description

ZEBRA Invento™ is a positioning chair designed for children and adolescents who need support in the process of sitting independently or in adopting and maintaining a sitting position.

The device combines the function of a seat and a lying position, which combined with the smooth adjustment of tilt-in-space function, supports the rehabilitation process. Adjustable headrest, footrests, chest and pelvic supports help maintain a correct lying, semi-lying or sitting position.

The functionality of **ZEBRA Invento™** makes the positioning chair meet all the expectations of even the most demanding users, significantly increasing everyday comfort, convenience and safety.

Size 1 available in december 2023 – ask the sales department for details.

Indications

- Cerebral palsy – various types and forms.
- Brain and spinal cord injuries with lower limbs paralysis and weakness, balance and motor coordination disorders with different etiology.
- Demyelinating disease (e.g. multiple sclerosis – SM).
- Meningocele.
- Muscles diseases with weakness and paralysis (dystrophy, myasthenia gravis, myositis).
- Genetic disorders.
- Degeneration and neurological diseases.
- Other condition connected with muscles weakness, paralysis and musculoskeletal disorders.
- After head and spinal cord injuries and spine trauma with paralysis and paresis.



Special features



Stable structure

The device ensures adequate stability and safety of user.



Eco-friendly materials

We use only human and environment friendly materials.



Quiet operation

Using of gas spring in the product ensures quiet operation.



Spacer upholstery

The device is equipped with breathable, hypoallergenic 3D spacer fabric. Now you can be sure that the skin stays dry for longer!

24⁺

24-month warranty

We offer 24-month warranty on the product. We also provide after sales service.



Easy to use

Intuitive way of adjustment and easy assembly.

2^w1

2 in 1

The product has 2 utility functions for use in 1 single device.



Tilt-In-Space

It is possible to change the tilt angle of the seat. Adjustable with a foot and a gas spring, thanks to which the caregiver has free hands during adjustment.



Multi-plane foot stabilization

Precise adjustment of the lower limbs is possible thanks to modern, independently adjustable foot stabilizers that allow for abduction, rotation, and depth adjustment (front-back).

Colors


Breathable upholstery - **ATMO™**



GRASS
Size 1*



RED
Size 2

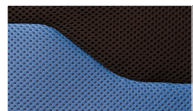

SEA
Optional

Size 1 available in december 2023 – ask the sales department for details.



Tilt-in-space adjustment

The gas spring allows for a smooth change of the seat inclination angle from (+20°) to (-25°). The adjustable foot leaves the caregiver's hands free, allowing for full control over the seat.



Lying position

Lying position – a unique function of the seat ensuring full comfort. With an extended backrest and/or deeper seat, additional filling pillows guarantee comfortable back support.



Additional filler cushions

Special cushions that ensure comfortable back support when extending the backrest and using the recline function


Standard equipment



+ Additional filler cushions



Important measurements

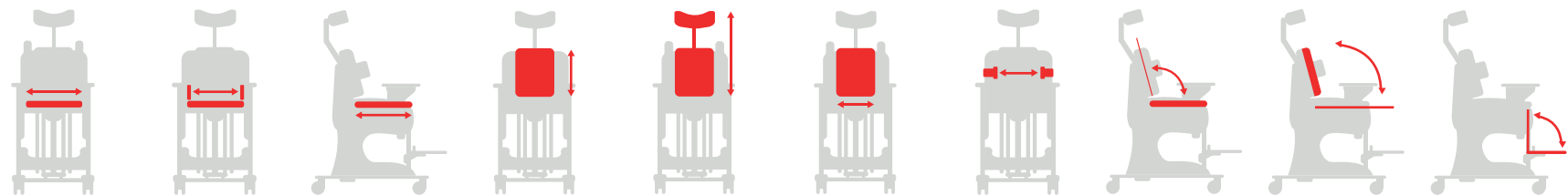

Sizes	Seat width	Seat width – with supports	Seat depth	Backrest height	Backrest height (without headrest)	Backrest width	Backrest width (with supports)	Seat tilt angle	Backrest tilt angle	Footrests tilt angle
Size 1* *	**Size 1 available in december 2023 – ask the sales department for details.									
Size 2	38 cm	18 - 26 cm	18 - 30 cm	32 - 42 cm	42 - 67 cm	40 cm	18 - 30 cm	-20° - +25°	90° - 180°	90° - 180°



Size	Indicative user height*	Maximum user weight
Size 2	80 - 130 cm	45 kg

* This measurement cannot be the main parameter while choosing the product size.
Note: Measurement tolerance range is +/- 1 cm / ½ in.
Note: Approximate values, depending on the product configuration.

Weight and dimensions

Size	Width	Length	Height	Weight
Size 2	56 cm	100 cm	95 cm	20 kg



Scan the QR code,
to see all the accessories.

Akcessories




Tray with edge

It allows you to play and perform therapy exercises. Perfect for working or taking meals.



Soft cover for tray

Soft cover protects the user from hitting the tray. Indicated for patients with impaired limbs coordination.

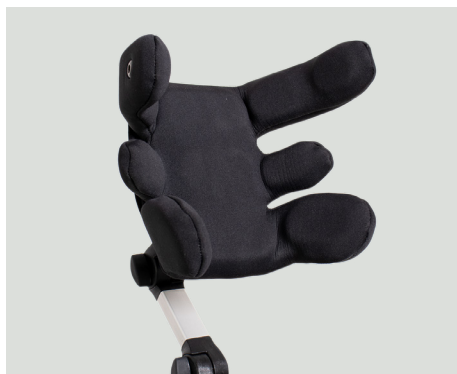

VisoPro™ headrest

The design of the headrest allows you to control head positioning at three levels. The individual parts of the headrest help protect the head from unwanted rotation and from falling to the sides.



Dynamic head supporting belt

Attached to the headrest using an elastic, dynamic tape with a velcro fastener. Used by people with limited head control to prevent it from drooping. Recommended for users with increased muscle tension. Special silicone elements prevent the belt from slipping.



Lower limbs supporting cushions

Mounted to the footrests. They provide additional support for the lower legs. Independent height adjustment allows them to be adjusted also to people with asymmetry of the length of the lower limbs.



Scan the QR code, to see all the accessories.